



Bedside Call Unit

ZKRUNCIPWL10100

ZKR

The Bedside Call Unit is used in patient rooms. There are backlit call and cancel buttons on the unit. In an emergency, a patient uses the call button to make an emergency call which appears as an alert on the Nurse Control Panel. Typical locations for this unit are on the walls of patient rooms and living areas.

- Handsets work through connection to this unit via RJ45 port on the front
- Can be flush mounted and has an ergonomic design
- RJ45 socket on the back provides easy installation
- Made of high quality PC-ABS-FR Raw Materials
- Call and Reset buttons are backlit and easy to notice in the dark

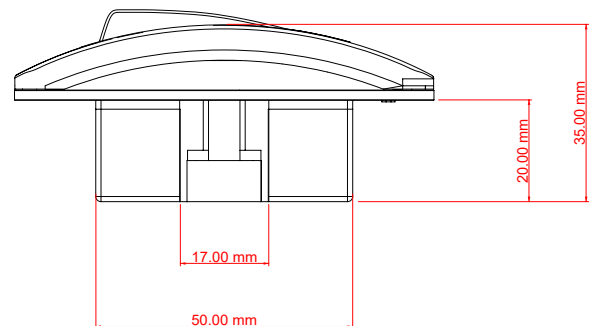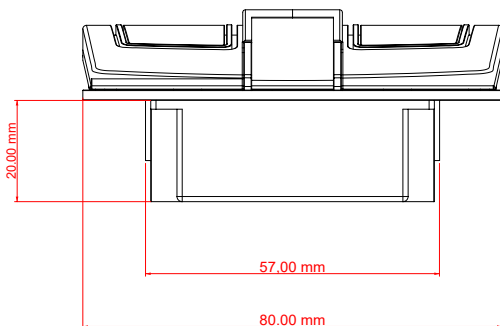

Technical Specification

ZKRUNCIPWL10100

<i>Weight</i>	56 g
<i>Width</i>	80 mm
<i>Depth</i>	35 mm
<i>Height</i>	80 mm
<i>Mounting Dimensions</i>	W-55mm D-20mm H-55mm
<i>Processor</i>	ARM
<i>Operating Voltage</i>	12V
<i>Operating Temperature / Humidity</i>	-10°C +60°C / % 30 %80 RH
<i>Raw Material</i>	PC-ABS-FR
<i>Buzzer Speaker</i>	2400Hz 65 dB
<i>Handset Connection</i>	1xRJ45 in front
<i>Ports and Slots</i>	1xRJ45 for connection to Room Control Unit
<i>Color</i>	White
<i>Back Light</i>	LED